

Sani-Cloth® Chlor: Training Poster

Two in one Cleaning & Disinfectant wipes

For non-invasive medical devices; non-porous hard surfaces.
It can be used to pre-clean or decontaminate medical devices prior to sterilisation. Latex Free. For professional use only.




Preparation



+ Remove lid, peel away and dispose of the film.



+ Take the first wipe and thread it through the opening in the canister lid.

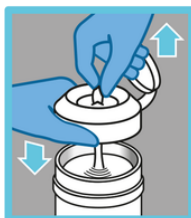

+ Replace the lid back onto the canister.



+ Dispense the wipe and open.

Cleaning & Disinfection



- + Wipe surface in an S-motion from top to bottom.
- + Ensure the surface remains visibly wet for the required contact time.
- + Use an additional wipe if needed. Allow surface to air dry.
- + Use one wipe per surface. Do not reuse wipes as this could impact the effectiveness or lead to cross-contamination.

Disposal



- + Dispose of wipe after a single use.
- + Packaging is recyclable.



Other information:

- + Coverage approx. 1 wipe per 1m²
- + Do not use if wipe is dry
- + For Professional Use only
- + Wear protective gloves

Compatibility



- + Compatible with over 900 medical devices and the most common materials in healthcare facilities.