

The Chest Seal you can trust



Versatile design for extreme situations

- Compact, fits IFAK without folding, to avoid creases.
- Strong adhesion in extreme conditions.
- Two per pack for ballistic wounds.
- Simple no-valve design.
- Vent with large, easy to operate tab.

From the makers of **CEL+X™**

actual size


Contact us at www.FoxSeal.com



Foxseal™

Foxseal™ Chest Seal is a pioneering occlusive dressing for the treatment of open chest wounds as a tool in prevention of pneumothorax. Other dressing options are often folded when carried in kits, leading to a permanent crease which can prevent a true seal. Foxseal is a compact size to fit neatly into pouches, to avoid creasing. Two seals per pack allow entry and exit wounds to be treated. Patient monitoring is vital in these injury types and the design is a simple valve-free option, for simple use. The chest can be vented with the easily-gripped tab and the adhesive re-seals.

Foxseal provides dependable adhesion in the most challenging situations.

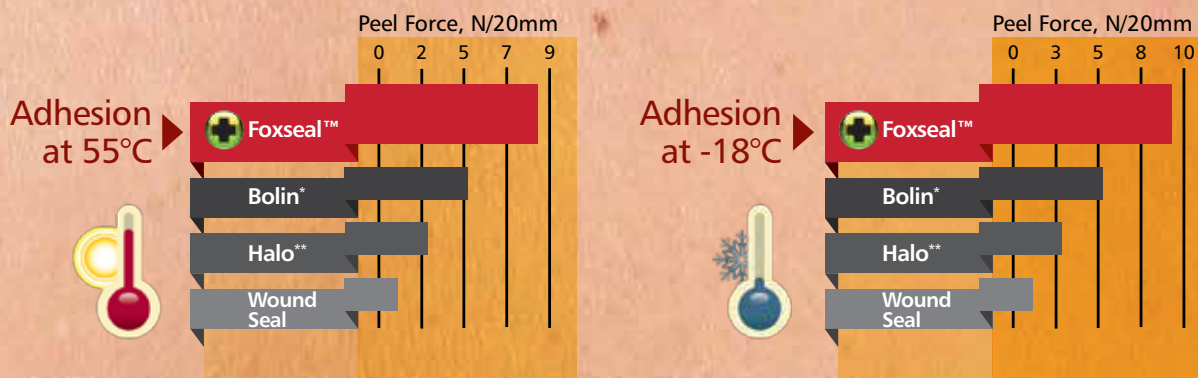

Foxseal™ information

- Sterile.
- Latex-free.
- Two seals per pack.
- FDA and CE listed.
- 4 tear notches for opening.
- Injuries up to 15cm maximum length.
- Adhesion performance maintained at -18 and +55°C.
- Compact pack fits IFAK without folding.
- Valve-free.
- Corrugated, high visibility peel tab for venting.
- Re-sealing adhesive.

Foxseal Chest seal is a latex-free single use dressing. Supplied sterile. Two chest seals per pack.

Ordering Information

Foxseal™:
One pouch containing two
Foxseal Dressings: FG08814441

Medtrade Products Ltd.

Electra House, Crewe Business Park, Crewe, CW1 6GL UK
T: +44 (0)1270 500019 F: +44 (0)1270 500045

Contact us at www.FoxSeal.com

MT/14/008

