



MOTION STRETCHERS

2 & 4 Section

Advanced ergonomic transport stretchers

Transport Stretcher

MOTION & MOTION PLUS



Structure

Epoxy coated steel tube frame.

Two or four section platform made of damp-proof and impact-resistant HPL material, resistant to humidity and shock.

Motion Plus (4 section) hydraulically operated.

Comfort

High-density polyurethane construction, ergonomically designed for optimum comfort.

Hydraulically operated backrest and height adjustment from 66–97.5 cm.

Hydraulic Trendelenburg and reverse Trendelenburg positions (+/-18°).

Foot pedals at either side of the stretcher control height, Trendelenburg and reverse Trendelenburg adjustment.

Easy Cleaning

Mattress upholstered with antimicrobial, anti-static covering. Easy to clean and maintain.

Smooth surfaces and removable panels allow easy, fast cleaning and disinfection.

Ease Of Handling

Ø 200 mm castors with single-point braking system, controlled by foot pedals on each corner of the stretcher.

The fifth wheel improves manoeuvrability and control; allows 360° rotation on its axis.

Ergonomically designed folding handles for easier greater control when manoeuvring.

Safety

X-ray translucent platform 10cm height mattress complies with BS 7176 (Medium Hazard).

Integrated siderails help secure patient. Their lowest position does not exceed the height of the mattress, making transfers simpler and safer.

Buffers located on each corner of the stretcher.

ABS base cover protects the stretcher mechanism and facilitates easy cleaning.

Conforms to Regulation (EU) 2017/745.

Specifications


MOTION & MOTION PLUS

2 Section stretcher: MOTION

Overall dimensions	218 x 82cm
Patient surface	190 x 72cm
Height adjustment	59–88cm
Back section dimensions	62 x 71cm
Leg section dimensions	66 x 116cm
Backrest tilting	90°
Weight w/o accessories	100kg
Safe working load	300kg
Trendelenburg / Reverse	±18°

4 Section stretcher: MOTION PLUS

Overall dimensions	210 x 82cm
Patient surface	190 x 72cm
Height adjustment	59–88cm
Back section dimensions	62 x 71cm
Leg section dimensions	62 x 34cm
Foot section dimensions	62 x 47cm
Backrest tilting	90°
Legrest tilting	35°
Weight w/o accessories	100kg
Safe working load	300kg
Trendelenburg / Reverse	±18°

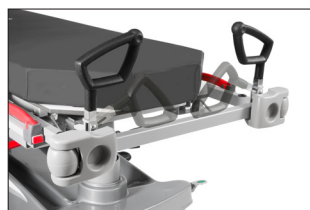

www.medi-plinth.co.uk

7–11 Holywells Road, Ipswich, Suffolk IP3 0DL – Call: +44 (0) 1473 212010