

**Baxter**



# Bravo Mini

---

24-Hour Ambulatory  
Blood Pressure Monitoring System

**SUNTECH**

# The SunTech Bravo Mini

The **Bravo Mini** system from SunTech Medical is a 24-hour ambulatory blood pressure monitoring (ABPM) device.



Data obtained from ambulatory blood pressure monitors is designed to be accurate and useful for managing a wide variety of hypertensive situations including white-coat hypertension.

## Clinically Validated

The **Bravo Mini** system is independently tested to meet the accuracy and performance requirements of the British Hypertension Society, the European Society of Hypertension International and the ISO 81060-2 standards.

## SunTech Motion Tolerance

The **Bravo Mini** device utilizes SunTech proprietary motion tolerant technology, to obtain reliable readings during motion.

## Smart Inflate

Smart inflate algorithm that bases measurement inflation on previous measurements to help protect against over-inflation and provide patient comfort.



## Orbit Blood Pressure Cuffs

SunTech's patented **Orbit** cuff is specifically engineered for ambulatory blood pressure monitoring applications. This cuff is designed with a soft sleeve that gently hugs the arm to help maintain cuff placement and provide patient comfort for 24 hours.



## AccuWin 4 SE Software

**AccuWin 4 SE** is a simple, **Windows** software based program that allows for configuration, analysis, interpretation and reporting of ABPM studies.



### VISUAL BLOOD PRESSURE RESULTS

The unique **AccuWin** report makes communication with patients simple with red-green graphical representation of blood pressure study data.


### AUTOMATIC INTERPRETIVE SUMMARY

Decipher data, creating interpretive summaries, according to AHA, ESH and Paediatric guidelines, for:

- Automated analysis of ABP levels for normality
- White-coat syndrome
- Overnight dipper status

## Client – Server Software

Access patient information from across your network with patient data securely stored in a networked server.

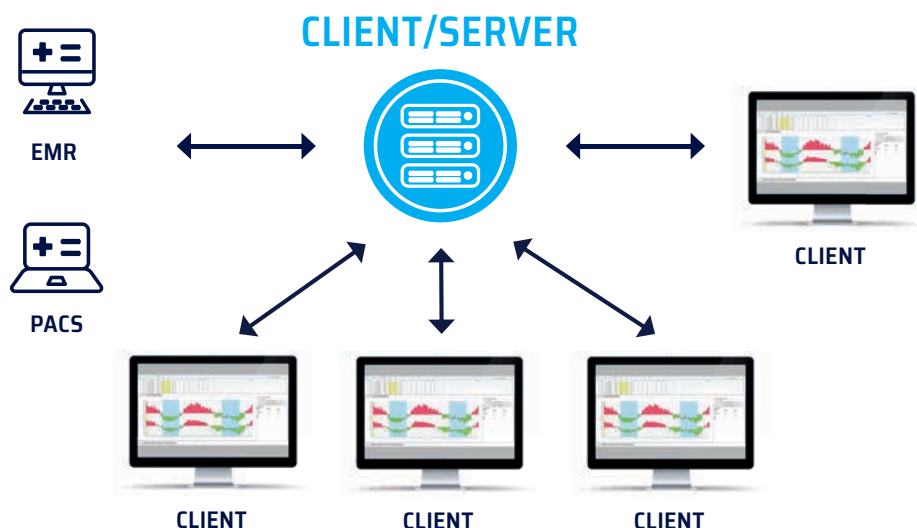

### SELECT CONFIGURATION AND ACCESSORIES

|        |                                                                                                                                            |
|--------|--------------------------------------------------------------------------------------------------------------------------------------------|
| 109700 | <b>Bravo Mini</b> ABPM Kit with software, cable, pouch, shoulder strap, belt, batteries, and 2 BP Cuffs (cuff sizes 26-34 cm and 32-44 cm) |
| 422850 | <b>AccuWin 4 SE</b> Software, USB drive                                                                                                    |
| 422851 | ABPM USB-C Cable                                                                                                                           |
| 422853 | #1 <b>Orbit</b> Comfort Cuff (18-26 cm)                                                                                                    |
| 422854 | #2 <b>Orbit</b> Comfort Cuff (26-34 cm)                                                                                                    |
| 422855 | #3 <b>Orbit</b> Comfort Cuff (32-44 cm)                                                                                                    |
| 422856 | #4 <b>Orbit</b> Comfort Cuff (42-55 cm)                                                                                                    |
| 422857 | ABPM Pouch                                                                                                                                 |
| 422859 | ABPM Shoulder Strap                                                                                                                        |
| 422860 | ABPM Belt                                                                                                                                  |

## Learn More About the Bravo Mini System, Contact Your Baxter Sales Representative.

Manufacturer: SunTech Medical, Inc. 5827 S. Miami Boulevard, Suite 100 Morrisville, North Carolina 27560-8394 USA. EU Authorised Representative: Emergo Europe, Westervoortsedijk 60, 6827 AT Arnhem, The Netherlands. Medical device (Regulation 2017/745): Class IIa. Notified Body Intertek 2862.

The Bravo Mini System is a non-invasive oscillometric ambulatory blood pressure monitor that is intended to be used with AccuWin 4 SE, a PC-based computer program for the recording and displaying of up to 250 measurements of systolic and diastolic blood pressure and heart rate. It is intended for use as an aid or adjunct to diagnosis and treatment when it is necessary to measure an adult and paediatric (> 3yrs.) patient's systolic and diastolic blood pressures over an extended period of time. The system is only for measurement, recording, and display to assist a licensed physician in making a diagnosis.

For safe and proper use of the product mentioned herein, please refer to the appropriate Operator's Manual or Instructions for Use. This medical device is a regulated health product which, pursuant to such regulation, bears a CE mark. Baxter recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical device. The personnel of healthcare establishments are responsible for the proper use and maintenance of this medical device. Baxter reserves the right to make changes without notice in design, specifications and models.

### Hillrom.Baxter.eu/en-GB

Baxter Healthcare Ltd  
Spur A, 2nd Floor, Reading International Business Park, Basingstoke Road, Reading, RG2 6DH

Baxter is a registered trademark of Baxter International Inc.  
AccuWin, Bravo, Bravo Mini, Orbit, SunTech and SunTech Design are trademarks of SunTech Medical, Inc.  
Windows is a trademark of Microsoft Corporation.

UKI-FLC220-260001 (v1.0) 07/2026