



TAPERED VACUUM MATTRESS CDK




FEATURES

- **Non absorbent and self-extinguishing** robust fabric
- **Tapered shape** to improve adherence to the patient
- Fitted with a **reinforced perimeter rope to attach the patient fixation belts**
- **9 rubber sleeves** improve patient transport operations
- The **top side is orange** and the **bottom side is grey**
- **Extremely resistant on/off valve** made of chromed brass with a built-in fitting for rubber tubes
- CDK allows for preliminary **X-ray diagnostics**.

CDK CDK is a robust and tapered vacuum mattress that features 9 soft carrying handles. It is designed to offer an optimal immobilization of the trauma patient also during transport.

The lower part of CDK is fitted with a rope reinforced with Santoprene® handgrips to ensure an easy handling system. On the top part the combination of the rope and the handles serves as a strong loop to attach the patient fixation belts. The top side is orange and the bottom side is grey.

The on/off valve made of chromed brass with a built-in fitting for rubber tube is resistant to high and low temperatures. Radio transparency allows for preliminary X-ray diagnostics.

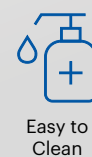

TECHNICAL DATA

| | |
|-----------------------------|------------------------------------|
| Length | 215 cm |
| Width | 90 cm |
| Weight (not including pump) | 4,1 kg |
| Handles | 9 |
| Load capacity | 150 kg |
| Materials | PP, PVC, Nylon, Brass, Polystyrene |
| Colour | Orange/Grey |

ACCESSORIES



STO0593A

K BELT 3

Set of 3 belts STX 597 with bag

- Materials: Nylon
- Plastic buckle
- Length min/max 40/190 cm
- Belt width: 5 cm
- Kit weight: 450 g



QM22120A

QMX120

Aluminium pump for vacuum mattresses

- Materials: Al, PVC, Steel
- Dimension: 50x14x8 cm
- Weight: 1,14 kg
- Tube length: 100 cm



QM22109A

QMX109

Pump in ABS for vacuum mattresses

- Materials: ABS, PE
- Dimension: 21,5x10x44,5 cm
- Weight: 0,6 kg
- Tube length: 100 cm



QM22170A

QMX 170

Bag for transport and storage of vacuum mattresses

- Materials: PVC, Nylon
- Dimensions: 103 x 87 cm
- Weight: 1,10 kg
- Fitted with zip closure



QM22199A

QMX 199

Repair kit for vacuum mattresses or vacuum splints

- 4 PVC patches
- 1 tube of glue (5 g)
- One patch size: Ø 5 cm
- Overall weight 14 g

Rev.0 (28/08/2023)

©Copyright - This document is property of SPENCER ITALIA S.r.l and can't be reproduced, used and disclosed without the written approval of SPENCER ITALIA SRL.