


Solite Pro

Community Bed Range

The highest specification
product in its class with
7 years warranty

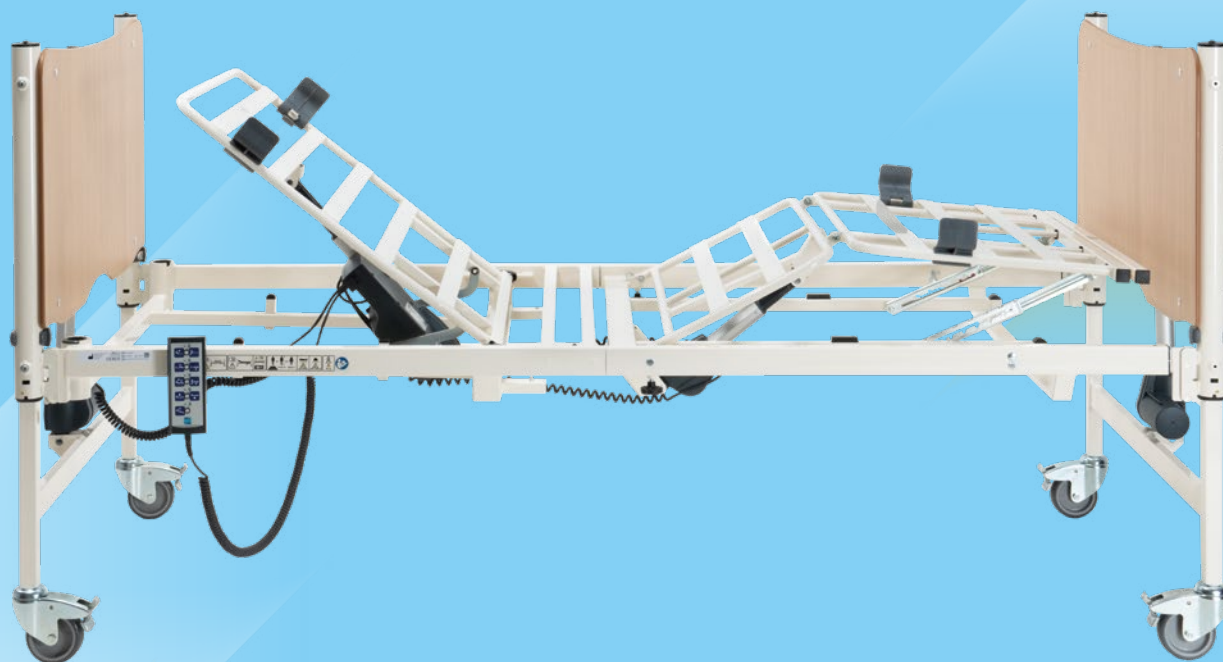

Maximise care, minimise risk

Constantly developed and improved over the years to provide the optimum combination of client comfort, versatility, functionality and performance.


The Solite Pro is a complete redesign of the Solite UK, now offering additional functionality whilst retaining an attractive domestic look and feel. This electrically operated, 4 section profiling bed features a built-in auto regression function, and transport stand storage locators. The optional patented side rails can be added for complete in-bed safety whilst ensuring easier ingress and egress for the patient.

Verified to IEC 60601-1 General Medical Electrical Equipment standard, IEC 60601-2-52 Medical Bed standard and IEC 60601-1-11 Medical Electrical Equipment standard specific to homecare requirements



Product Overview

Enhancing Lives Together

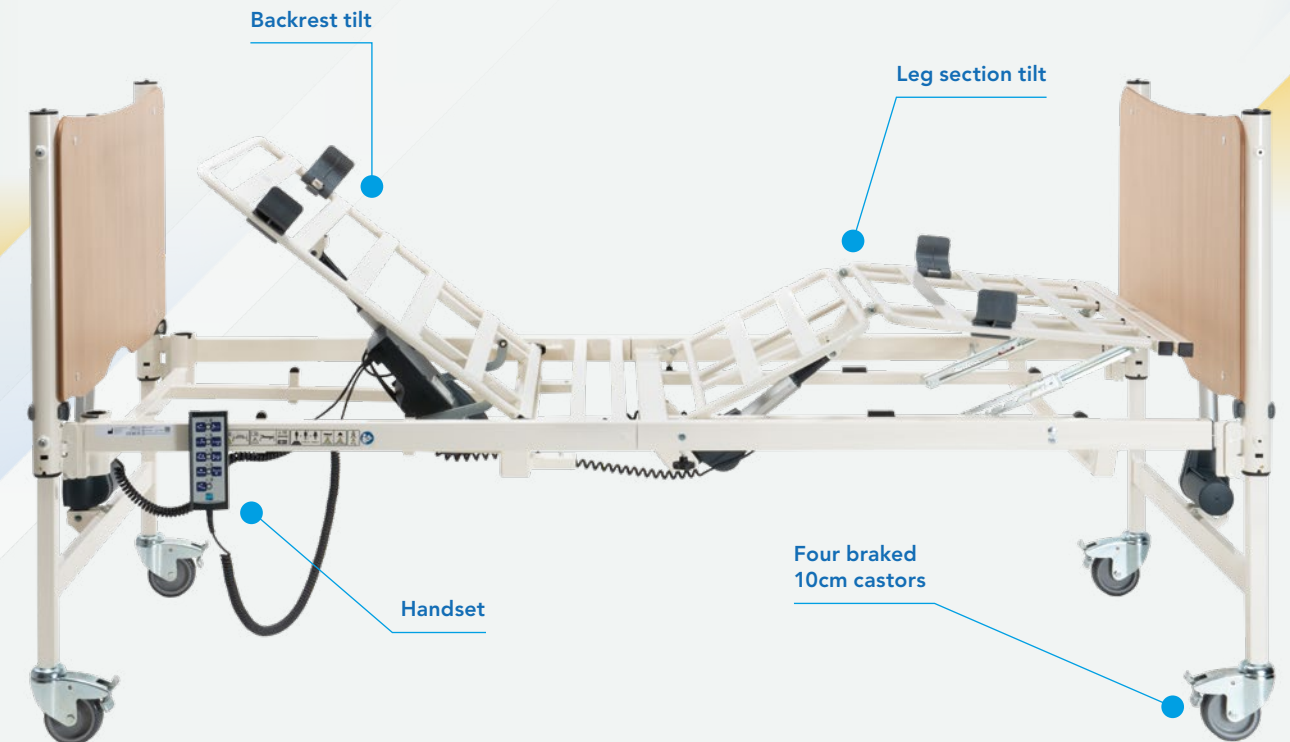

Key Features

Individual sections strengthened for extended performance and reliability

Improved auto regressing backrest prevents pinching in the sacral area

Four braked castors for manoeuvrability and safety

Easy-use handset controls all functions – includes lockout function

Tool Free Assembly

Single Driver installation

Knee raise function

12° foot down tilt with head down tilt available as an optional extra

In-built transport stand storage locators below the mattress platform



Electrically profiling bed with back rest and knee break function, variable height control and auto contour



Slatted Bed Platform – Ensures easy cleaning, durability, and improved infection control



Built-in integral leg extension for taller patients, including auto stop for added safety

Solite Pro Community Bed Range



Solite Pro
Engineered for comfort, built for reliability



Solite Pro Low
Innovative design, effortless care




Solite Pro UltraLow
Ultra-low design, ultimate fall protection

Technical Details & Warranty

drivedevilbiss.co.uk


Mattress Platform Height Range

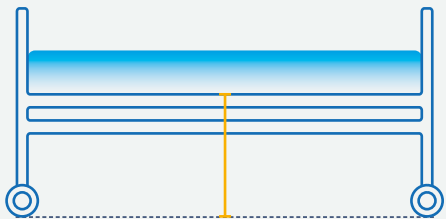

39.5 - 81.5cm
Solite Pro



29.5 - 71.5cm
Solite Pro Low

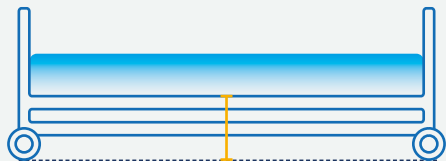

7.5 - 80.5cm
Solite Pro UltraLow

Advanced Profiling Functions



Auto Contour
A pre-programmed profiling function that raises and lowers the backrest and knee break simultaneously. This helps maintain alignment and reduces sliding down the bed.



Cardiac Chair Position
Creates a semi-seated, chair-like posture within the bed, combining the Auto Contour with optional Reverse Trendelenburg tilt. Ideal for patients who remain in bed for extended periods, it encourages interaction, self-feeding, all while maintaining comfort and support.



Auto Regression
As the backrest rises (0°-70°), the base slides slightly backward to increase space in the pelvic area. This reduces shear, enhances comfort, and mimics the spine's natural movement and is especially beneficial during backrest elevation.



Reverse Trendelenburg
Provides a foot-down tilt of up to 12°, elevating the head end of the bed. This position supports improved breathing, digestion and social engagement, while offering an alternative lying posture to reduce discomfort or pressure.



7 Years of Warranty
When you choose our products, you're not just investing in exceptional British craftsmanship, backed by over 130 years of UK manufacturing expertise, you're also securing peace of mind for years to come. This product comes with a guaranteed 7-year warranty, which can be extended to 10 years when a Sidhil service contract remains in place.



Why Choose
Drive DeVilbiss?

Wherever care takes place, we're here to help you deliver it with confidence.

At Drive DeVilbiss, we understand that care isn't just about clinical need. It's about supporting daily routines, restoring confidence, and helping people feel safe, comfortable and in control.

We've supported hospitals, care homes and community services across the UK, but what drives us is something more personal. People. The everyday moments that matter. Whether it's enabling a better night's sleep or helping carers move safely, we focus on the small things that improve lives in a big way. That's why everything we design is built with care, with clarity, and with the belief that everyone deserves to be supported.



Enhancing Lives Together



Here for you, together...

As part of a global healthcare group rooted in the everyday realities of UK care, we work with providers and professionals who share our values. These are people who listen, understand, and are committed to delivering the right support at the right time.

When you choose Drive DeVilbiss, you're backed by experience, by people who care, and by a network that's here to support you now and in the years ahead.

Trusted expertise, close at hand

We are trusted by care professionals, care home groups, and the NHS. That same clinical insight goes into every product we design and deliver. Our dedicated UK-wide team has the knowledge and experience to guide you through product selection, installation, and ongoing support. Wherever you are based, help is never far away.

130 Years of UK Manufacturing

Our Halifax production site directly employs 128 people and supports a network of UK suppliers, sustaining an additional 75 jobs. In total, our UK manufacturing contributes £1.52 million in tax revenue every year. This is enough to fund the equivalent of 2,916 NHS hospital beds. By sourcing 70% of materials locally, we invest in jobs, skills, and the future of British industry.

Better for the planet

Choosing UK manufacturing significantly reduces our environmental impact. Transporting 20 tonnes of medical beds from China to the UK creates over 28,000kg of CO2 emissions. Producing and delivering the same quantity from Halifax generates just 552kg of CO2, which is a 98% reduction. By avoiding long-distance shipping, we also cut other harmful emissions such as nitrogen oxide and particulate matter.

Codes and Specifications

Solite Pro

Width: 96.5cm (101.5 with side rail channels)
Height: 39.5 - 81.5cm (to mattress platform)
Length: 223cm (239.5 when fully extended)
Weight: 85kg
Max. user weight: 191kg (30 stone)

Solite Pro Low

Width: 96.5cm (101.5 with side rail channels)
Height: 29.5 - 71.5cm (to mattress platform)
Length: 223cm (239.5 when fully extended)
Weight: 85kg
Max. user weight: 191kg (30 stone)

Solite Pro UltraLow

Width: 95.1cm
Height: 7.5 - 80.3cm (to mattress platform)
Length: 239.7cm
Weight: 115kg
Max. User Weight: 160kg (25st)



Scan for full features and resources.

Drive DeVilbiss Healthcare Ltd. Sidhil Business Park,
Holmfild, Halifax, HX2 9TN. Item Code: LL/SOLITE/PRO/8pp/L

We have a continuous product improvement policy and consequently reserve the right to amend design and specification without prior notice. All sizes and weights stated are nominal. All errors and omissions excepted. Due to the limitations of photography and the printing process colours shown may not be 100% accurate. *Please note: To fit this accessory then side rail channels must be present. Accessories should be fitted by a competent technician



drivedevilbiss.co.uk

Solite Pro Range

Solite Pro Bed: SOLITE/PRO
Solite Pro Low Bed: SOLITE/PRO/LOW
Solite Pro Low with Battery Back Up: SOLITE/PRO/LOW/BB
Solite Pro Ultra Low Bed: SOLITE/PRO/ULOW
Side Rail Channels: SOLITE/PRO/ADD/CHANNELS
Side Rail Channels: SOLITE/PRO/
Factory fitted to bed: CHANNEL/FITTING

Side Rails, Pads &, End Pads

Full Length Side Rails*: SOLITE/PRO/SR/A2
Extended Length Side Rails*: SOLITE/PRO/SR/EXT
Side Rail Height Extensions: SOLITE/PRO/SR/HE
Grange Side Rails: GRANGE/SR
Standard Side Rail Pads: SOLITE/PRO/SR/PAD
Side Rail Length Extension Pads: SOLITE/PRO/SR/EXT/PAD
Side Rail Height Extension Pads: BED01/SR/HEXT/PAD
Standard Side Rail Mesh Cradle: SOLITE/PRO/EPS
Grange Side Rail Pads: 1316/PS/CREAM
Solite Pro End Pads: SOLITE/PRO/HF/PAD

Solite Pro Accessories:

Transfusion Pole: BED/DP
Lifting Pole: SOLITE/LASER/LP/1
Cromer Cantilever Bed Cradle: 1369
Spare Lockable Handset: CONTROL/127
Lockable Handset, Head Down Tilt Function: CONTROL/038
Solite Pro Side Rail Transport Accessory: SOLITE/PRO/SR/TS
Handset Locking Key: CONTROL/049

Solite Pro Ultra Low Accessories:

Bed Extension (+16.5cm): ULOW/BED/EX
Mattress Extension Infill: ULOW/MAT/EX/INFILL
Raised Height Transport Stand: ULOW/HE/TS
Pump Bracket: ULOW/PUMPBKT
10 Button Handset: CONTROL/130
Dedicated Mattress: ULOW/MAT
Foldable Crash Mat 150mm: ULOW/MAT/FALL
Mesh Safe Side Rails: ULOW/MESH/RAIL

Grab Handles:

Elland Standard Grab Handle: 1210/GRAB/A2
Elland Short Grab Handle: 1216/GRAB/A2
Ultra Low Grab Rail: ULOW/GRAB/RAIL