



Temperature controlled AED cabinet ...Any place ...Any environment



Lockable keypad



Unlockable version



Personalised brand area



Temperature controlled



Intelligent surround lighting



LIFEVAULT™

Designed to keep **your** AED rescue ready.

Automatic heating and cooling system maintains a safe operating temperature whatever the weather and environment.

Too Hot

Heat can also drain batteries quickly. Electronics can fail as can components such as the pads which can dry out and fail to make a strong contact with the skin.




Too Cold

If the unit is too cold batteries can drain quickly, pads will not adhere, and the AED can malfunction rendering it useless in an emergency.



LIFEVAULT™

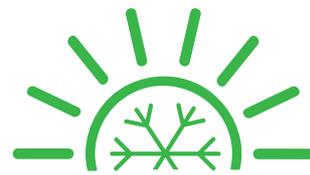

TEMPERATURE CONTROLLED AED cabinet

Easily accessed in an emergency

Available in keypad lockable and unlocked versions.

Why is temperature important?

Extremes in temperature can drastically affect the safe operation of the AED.

Designed for **every**

environment

-10°C to 50°C



Strong. Secure. Intelligent.

- Designed to fit ALL AEDs currently available.
- Automatically cools or heats to maintain a safe operating temperature between 10°C and 40°C. The cabinet automatically cools at an internal temperature of 35°C and heats up at 12°C.
- Intelligent LED lighting system. Light sensor triggers interior white light in low light conditions and optional external green LED surround light.
- Strong and sturdy vandal proof construction.
- Weatherproof IP66 construction.
- Securely lockable with keypad - unlocked version also available.
- Double walled insulated cabinet for energy efficiency.

LIFEVAULT AED CABINET



7900 LIFEVAULT TEMPERATURE CONTROLLED AED CABINET Lockable door

SIZE: 56.2cmH x 45.3cmW x 27.9mD
UNIT: 1
MOQ: 1



7901 LIFEVAULT TEMPERATURE CONTROLLED AED CABINET Unlockable door

SIZE: 56.2cmH x 45.3cmW x 27.9mD
UNIT: 1
MOQ: 1

Suitable for the **MEDIANA** Automated External Defibrillators

MEDIANA defibrillators and accessories are available at www.medianadefib.co.uk

FREE E-Learning
course



MEDIANA A15

Automated External Defibrillator

- Adult / Child setting with one switch
- Suitable for workplaces, schools and clubs
- Automatic switch-on when lid opens



MEDIANA A16

Fully Automatic Defibrillator

- 3 pre-loaded languages - English / Welsh / Polish
- Fast shock delivery
- Rated to IP55



+44 (0) 08456 448808 | sales@reliancemedical.co.uk | www.reliancemedical.co.uk

RELIANCE MEDICAL
"Leading First Aid"



REV: 02.25