



THE ULTIMATE IN PASSIVE REWARMING TECHNOLOGY

Blizzard with Reflexcell™ technology traps warm air in the blanket's cells, providing life-saving insulation. The elastic inner bands keep the blanket close to the body, reducing convection, while the silver reflective surface blocks radiated heat loss. The outer layer is completely waterproof and windproof.


Blizzard's insulating material is rated in Togs, a measure of thermal resistance. Our 2- and 3-layer systems, with a 5.5–7.5 Tog rating, are equivalent to a medium-weight sleeping bag.

The Journal of Trauma: Injury, Infection, and Critical Care stated that Blizzard was more effective and outperformed competitors' active rewarming products in a 2010 study.¹ Blizzard with Reflexcell is the standard wherever lightweight, compact thermal performance is essential.

1 Preventing Hypothermia: Comparison of Current Devices Used by the US Army in an In Vitro Warmed Fluid Model. Allen, Paul B. DSc, PA-C; Salyer, Steven W. PA-C; Dubick, Michael A. PhD; Holcomb, John B. MD; Blackburn, Lorne H. MD. Author Information: The Journal of Trauma: Injury, Infection, and Critical Care 69(1):p S154-S161, July 2010. | DOI: 10.1097/TA.0b013e3181e45ba5.

HOW DOES Reflexcell™ TECHNOLOGY WORK?

- 1 Cells trap air
- 2 Elastic hugs body
- 3 Metal coating reflects heat

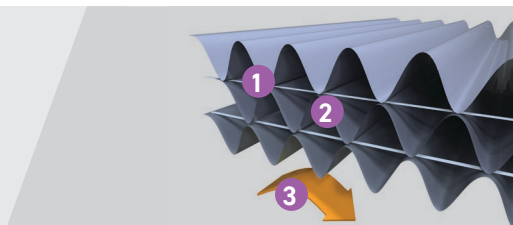

PASSIVE REWARMING FOR HYPOTHERMIA TREATMENT:

Blizzard 2-layer and 3-layer systems trap warm air to effectively treat hypothermia due to blood loss and environmental exposure.

STURDY DESIGN: Reinforced edges offer maximum strength, while a resealable self-adhesive closure allows for rapid patient access.

VERSATILE APPLICATIONS: Military, hospital, emergency service, search & rescue, outdoor enthusiasts, veterinary care, sports, and industrial applications.

PROVEN EFFECTIVE: CE and FDA approved.



visit safeguardmedical.com
call (US) +1 (855) 428-6074 / (Int'l) +44 (0) 1568 613942
email (US) info@safeguardmedical.com / (Int'l) sales.admin@safeguardmedical.com



Product name	PN	NSN	Color	Layers	Length	Packaged dims.	Product dims.
Blizzard Trauma Blanket	BPS-01	6532-01-524-6932	Green	3	Full	9.1" x 4.3" x 2" (23 x 11 x 5cm)	7.5ft x 6.5ft (2.3 x 2m)
	BPS-10	7210-99-958-4316	Orange				



Product name	PN	NSN	Color	Layers	Length	Packaged dims.	Product dims.
Blizzard IFAK Blanket	BPS-47	4240-99-334-1172	Yellow	2	Half	6.7" x 4.3" x 1" (17 x 11 x 2.5cm)	3.6ft x 5.1ft (1.1 x 1.56m)



Product name	PN	NSN	Color	Layers	Length	Packaged dims.	Product dims.
Blizzard Survival Bag	BPS-02	8465-99-433-4810	Green	3	Full	9.1" x 4.3" x 2" (23 x 11 x 5cm)	7.5ft x 6.5ft (2.3 x 2m)
	BPS-11	8465-99-958-4315	Orange				



Product name	PN	NSN	Color	Layers	Length	Packaged dims.	Product dims.
Blizzard Survival Jacket	BPS-03	8415-99-574-0495	Green	3	Half	5.9" x 4.3" x 2" (15 x 11 x 5cm)	4 x 2.6ft (1.2 x 0.8m)
	BPS-12	8415-99-852-5983	Orange				



Product name	PN	NSN	Color	Layers	Length	Packaged dims.	Product dims.
Blizzard Baby Wrap	BPS-23	--	Orange	3	Full	3.9" x 3.2" x 0.8" (10 x 8 x 2cm)	25.6" x 9.4" (65 x 24cm)



Product name	PN	NSN	Color	Layers	Length	Packaged dims.	Product dims.
Blizzard Adult Beanie Hat	BPS-24	--	Orange	3	--	4.7" x 4.3" x 0.4" (12 x 11 x 1cm)	7.9" x 10.2" (20 x 26cm)

