


 **Medisave**  
www.medisave.co.uk

**seca** authorised dealer

go to [www.medisave.co.uk](http://www.medisave.co.uk) for the full  
seca range

# seca 877

Flat scales for  
mobile use 

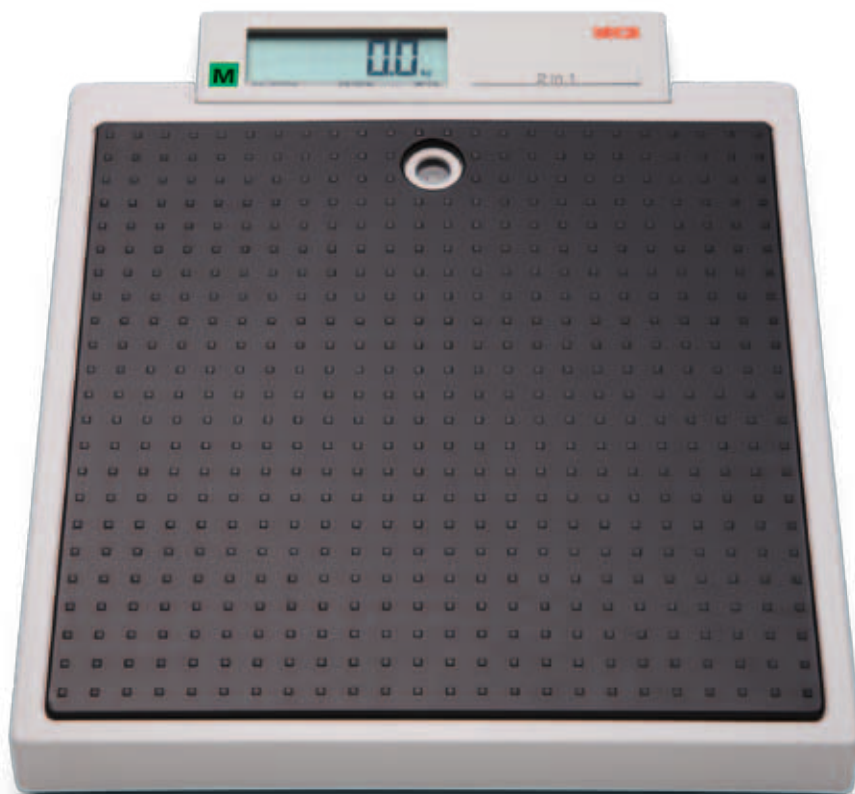

- Mother-Child function for trouble-free weighing of small children.
- Large leveling base for high stability.
- Easy to transport – scale itself weighs only 4,2 kg.
- Battery power for sites without electrical outlets.

## seca 877:

### A lightweight for heavyweight use.

In the seca 877 you've found the right solution if you're looking for a lightweight scale (only 4.2 kg) that's stable, compact, easy to clean and equipped with a Mother-Child function.



#### Light but stable.

The base is made of a mixture of ABS plastics and Plexiglas. This combination guarantees a material with outstanding surface hardness and impact strength. Furthermore, it is easy to clean with any commercially available disinfectant. Additional steel elements in the base make sure that persons up to 200 kg can weigh themselves easily.

#### Light but impact resistant.

ABS plastics and Plexiglas are not only very stable but are also especially light, a quality that is important for those interested in mobile use of the seca 877. Just as important for use away from home is the large leveling base which allows the scale to be positioned without any trouble regardless of the surface underneath it. It's easy to step onto the scale, thanks to the flat construction and slip-proof standing surface.

#### Light, but with Mother-Child function.

The scale is easy to operate even with the display and operating elements integrated in the base. Just a light tap of the toe on the scale surface turns on the scale. Infants and small children can be easily weighed while being held by the mother. All she has to do is weigh herself separately first. Then the Mother-Child function ascertains the tare of the weight.




Large display for easy reading.



Large leveling base on the bottom of the scale for safe placement.



Tip-on function for weighing without turning on the scales.

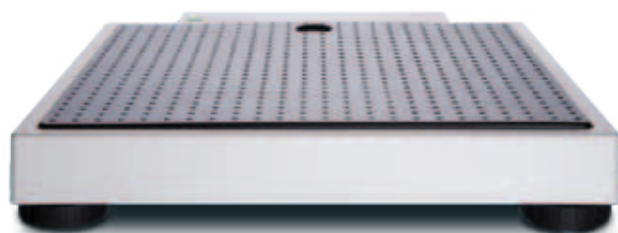

Flat construction and slip-proof surface for a comfortable foothold.

## seca 877

### Technical Data

- Weight capacity: 200 kg
- Scale division: 100g < 150 kg > 200 g
- Scale dimensions (WxHxD): 321 x 60 x 356 mm
- Weight of scales: about 4,2 kg
- Power: batteries
- Functions: 2 in 1 (Mother-Child function), Tip-on automatic for switching ON and OFF
- Optional: stadiometer seca 217 combined with adapter element seca 437, carrying case seca 414 (for combination with stadiometer and adapter element) and seca 421
- Calibration class: