

Nevershare 2mL

Luer lock syringes



EXCHANGE
SUPPLIES

Sterile 2 mL Luer lock syringe for single use. There are five different colour plungers: White; Blue; Green; Yellow; Pink.

Supplied in boxes of 100 syringes, which all have the same colour plunger. The boxes of 100 are supplied in outer cartons containing 3,200 syringes (see page 3 for details) or mixed cartons.

Product codes

S242, S243, S244, S245, S246, S247 (mixed)

Materials

Syringe body: Medical grade Polypropylene (PP 5250T)

Plunger gasket: Isoprene rubber (latex-free)

Syringe barrel and gasket lubricating oil:

Medical grade Silicone oil (DC360)

None of these materials contain recycled content

Sterilisation

Sterilised by: Ethylene oxide gas

Conforming to: ISO:11135 Sterilisation of medical devices
(validation and routine control of ethylene oxide sterilisation)

Time to expiry: 5 years from manufacture

Product conformance

Class IIa Medical Device

Fully compliant with ISO:7886-1: Sterile hypodermic syringes for single use

European Medical Device Directive 93/42/EEC (CE Class:1)

Manufacturing compliant with Quality Assurance Standards ISO:13485

Compatibility

All Luer fit needles

Manufacture and distribution


Nevershare syringes are manufactured exclusively for Exchange Supplies, UK.

1, Great Western Industrial Centre, Dorchester Dorset DT1 1RD. www.exchangesupplies.org

Manufactured for	Shanghai Kindly Enterprise Development Group Co., Ltd.
Exchange Supplies by:	No. 658 Gaochao Road, 201803 Shanghai, People's Republic of China.

EC representative:	Shanghai International Holding Corp. GmbH (Europe). Eiffestraße 80, 20537 Hamburg, Germany. T: 0049-40-2513175. F: 0049-40-255726
---------------------------	--

KDL UKCA representative:	KDL – Overseas Business Development Director. Unit 14 Cowley Mill Road, Uxbridge, Greater London UB8 2DB.
---------------------------------	--

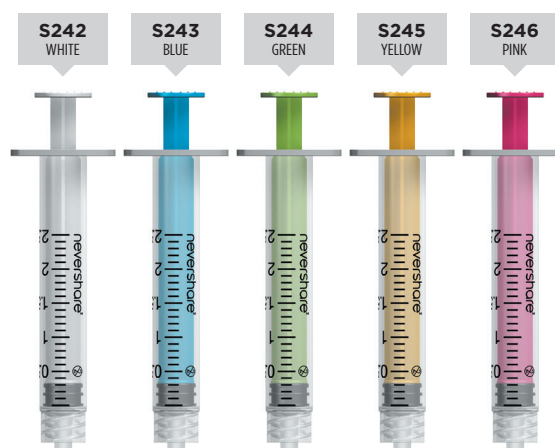

Graduations in
0.1mL increments

Packaging

Item	Weight (approx.)	Material(s)
Blister	4 g	Printed surface: Medical grade 63gsm coated paper Blister: 0.09m thickness medical complex film
Box 100	463g	E corrugated card
Carton 3,200	17.5kg	Uncoated 630gsm AB flute double wall card. Consisting of: 150gsm 3A Kraft Test Liner printed outer layer, with 120gsm and 100gsm A grade Senior corrugated and 140gsm B grade Kraft Test inner Liner. Recyclable.

None of these materials contain recycled content.



Blister
138.6mm x 39.5mm



Carton 3,200
Dimensions: 585mm x 440mm x 545mm



Box 100
210mm x 125mm x 140mm

S242

Blister GTIN: 06930197835519
Box GTIN: 16930197835516
Carton GTIN: 26930197835513

S243

Blister GTIN: 06930197835526
Box GTIN: 16930197835523
Carton GTIN: 26930197835520

S244

Blister GTIN: 06930197835533
Box GTIN: 16930197835530
Carton GTIN: 26930197835537

S245

Blister GTIN: 06930197835540
Box GTIN: 16930197835547
Carton GTIN: 26930197835544

S246

Blister GTIN: 06930197835557
Box GTIN: 16930197835554
Carton GTIN: 26930197835551

EXCHANGE
S U P P L I E S

1 Great Western Industrial Centre, Dorchester,
Dorset DT1 1RD | T: +44 (0) 1305 262244