

SAM
MEDICAL

SAM[®] Splint

SHAPE. SET. SECURE.

UNIVERSAL FLEXIBILITY


Touting unrivaled flexibility, SAM[®] Splint bends into any simple curve, becoming exponentially stronger and more supportive. SAM[®] Splint is globally acclaimed by emergency care providers, outdoor enthusiasts, the U.S. military, and even NASA.



FEATURED TECHNOLOGY

ADAPTABLE C-CURVE™ PRINCIPLE

SAM[®] Splint leverages the revolutionary C-Curve™ Principle to mold its foam and aluminum layers into structural arcs, strong enough for pre- or post-cast care, while remaining lighter and more transportable than traditional splints that rely on heavy, rigid materials. Bend it back into its original form to stow away or remold for a separate injury.

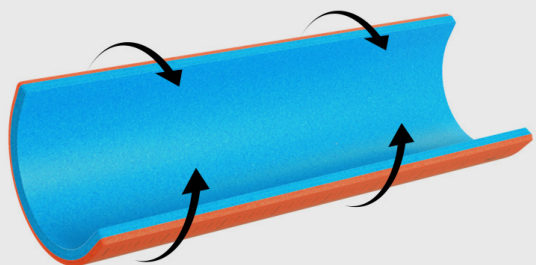

ENGINEERED FOR VERSATILE SUPPORT

1 CUT-TO-SIZE

SAM[®] Splint's two layers of closed-cell foam and its aluminum interior can be cut with ordinary scissors to adapt to any size requirements.¹

2 MOLDABLE ALUMINUM

Flat, O-temper aluminum inside the splint's outer layers molds easily, allowing for a wide spectrum of applications.

3 CLEANABLE EXTERIOR

Exterior closed-pore foam can be cleaned and disinfected with bleach and water for reuse.²

4 RADIOLUCENT MATERIAL


Medical professionals are able to take X-rays and CT-scans without removing the splint from the fractured or injured limb.

5 WEATHER RESISTANT

Designed to function across the spectrum of global temperatures – from the Sahara to Mt. Everest.

6 LIGHTWEIGHT DESIGN

Substantially lighter than plaster or fiberglass splints, SAM[®] Splint requires only a wrap or tape to secure an injured bone or stabilize a joint.

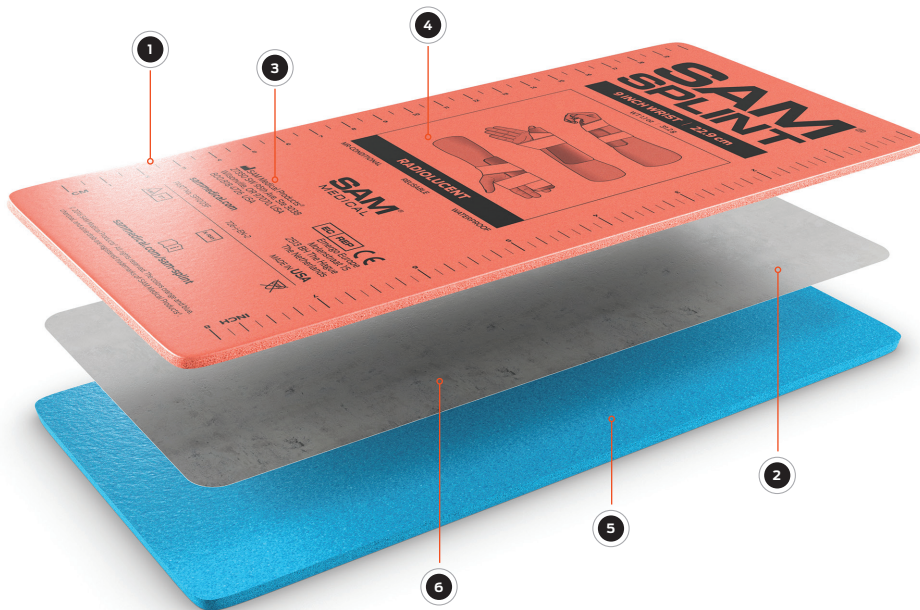

¹ Do not use a serrated blade. Ensure there are no sharp edges after cutting, and cover with tape.

² Care and maintenance guidelines are available at www.sammedical.com

MODELS



BLUE
ORANGE



CHARCOAL
CHARCOAL



OD GREEN
OD GREEN



A VARIETY OF
SIZES AVAILABLE

PRODUCT SPECIFICATIONS

BLUE / ORANGE MODELS*

36" Roll (SP506-OB-EN)
36" Fold (SP507-OB-EN)
36" XL (SP508-OB-EN)
18" (SP502-OB-EN)
9" (SP500-OB-EN)
Finger (SP510-OB-EN)
Industrial (SP504-OB-EN)

CHARCOAL MODELS*

9" (SP500-CH-EN)
18" (SP502-CH-EN)
36" Roll (SP506-CH-EN)
36" Fold (SP507-CH-EN)
36" XL (SP508-CH-EN)

OD GREEN MODELS*

36" Fold (SP507-OD-EN)

VET MODELS

36" VET (SP526-OB-EN-543)

GLOBAL SALES

U.S. Toll Free: 800.818.4726

U.S. Local: 503.639.5474

U.S. Fax: 503.639.5425

sammedical.com



Emergo Europe, Prinsessegracht 20,
2514 AP, The Hague, The Netherlands



SAM Medical Products[®]
12200 SW Tualatin Rd., Ste 200
Tualatin, OR 97062 USA



* Standard width 4.25 in wide, XL 5.50 in wide.

SAM[®] is a registered trademark in the United States and other countries.

©2019 SAM Medical Products. All rights reserved.