



TEST REPORT

Applicant: Guilin Refine Medical Instrument Co., LTD
Address: NO.8-3, Information Industrial Park, High-Tech Zone, Qixing District, Guilin, Guangxi, 541004, P.R.China

Report No: C80272046
Date: April 26, 2020

Attn: Chen Ju Dong
Sample Description: : Safety Goggle
Submitted sample(s) said to be: : 20 pieces
Item No.: : EP03
Device Type: : Safety Goggle
Material: : Polymer material
Frame Color: : Transparent
Lens Color: : Transparent
Lenses Type: : Uniform
P.O. No.: : Not Provided
Age Grade: : Not Provided
Supplier: : Guilin Refine Medical Instrument Co., LTD
Buyer: : Not Provided
Goods exported to: : Not Provided
Country of origin: : Not Provided
Date(s) of samples received: : April 9, 2020
Date of Test Period: : April 9, 2020 – April 26, 2020

REMARK

The Personal eye - protection test was performed at an UL approved subcontract lab.
 This report shall not be reproduced except in full, without the written approval of Shenzhen Precision Eyewear Testing & Inspection Services Co., Ltd.

Subcontract Report No.: PL 2004308

Notes: The results relate only to the items tested.

For and on behalf of :

UL VS Shanghai Limited Shenzhen Branch

Hedy Li



Hedy Li
 Consumer Engineering Manager
 Eyewear Testing Department

以下条款仅针对中国市场和社会:

——该报告仅向委托人发布、供委托人使用，不具有社会证明的作用。

This letter / report / certificate shall not be reproduced (except in full version) without the written approval of the UL Company. ("UL")

LETTERS / REPORTS / CERTIFICATES: Letters / Reports / Certificates of UL are issued for the exclusive use of the Customer to whom they are addressed. No quotation from reports / certificates or use of the UL's name is permitted except by UL's express written authorization. Letters / reports / certificates apply only to the specific materials, products or processes tested, examined or surveyed and are not necessarily indicative of the qualities of apparently identical or similar materials, products or processes. Reports / Certificates of UL do not relieve sellers / suppliers from their contractual responsibilities with regard to the quality / quantity of the goods in delivery nor do they prejudice the Customer's right to claim against sellers / suppliers for compensation for any apparent and/or hidden defects not detected during UL's random inspection or testing or audit. The liability of UL to the Customer in contract, tort (including negligence or breach of statutory duty) or howsoever, and whatever the cause thereof, (a) for any loss of profit, business, contracts, revenues, or anticipated savings; or (b) for any special indirect or consequential damage of any nature whatsoever, shall be limited to the amount of the fee paid in respect of the specific Work(s) which give rise to such claim. For Reports / Certificates, refer to the UL VS Terms and Conditions.

SZ-FAF-001 (2017-10-12)

UL VS Shanghai Limited Shenzhen Branch
 3-4F Qingyi Supermask Photoelectricity Building, No.8 Langshan 2nd Road,
 North High-Tech Industrial Park, Nanshan District, Shenzhen, P.R. China

优力检测质量检测(上海)有限公司深圳分公司
 广东省深圳市宝安区西乡街道铁岗社区铁岗山二路8号铁岗光电大厦三层。邮编
 电话(T): 755-26018990 / 传真(F): 755-26018626 / 网站(W): W: ul.com



Test Report

Report No.: PL2004308

Page 1 of 10

Applicant: Guilin Refine Medical Instrument Co., LTD

Address of Applicant: NO.8-3, Information Industrial Park, High-Tech Zone, Qixing District, Guilin, Guangxi, 541004, P.R.China

Date of Receiving Samples: Apr 13, 2020

Testing Period: Apr 14, 2020 to Apr 17, 2020

Description of Samples
The submitted sample and sample information was/were submitted and identified by/on behalf of client;

Sample Name: Safety Goggles
Model No.: EP03
Quantity: 12 Pairs
Material: Polymer Material
Scale No.: Not provided
Frame Color: Transparent
Lenses Color: Transparent
P.O. No.: Not provided
Supplier / Brand: Not provided
Buyer: Not provided
Goods exported to: Not provided
Country of Origin: China

Requested Standard:

EN 166: 2001	Personal eye - protection – Specifications
EN 167: 2001	Personal eye - protection – Optical test methods
EN 168: 2001	Personal eye - protection - Non-optical test methods

Results/Remarks: Please refer to the following page(s).

Issued by stamp

Date of Issued: Apr 24, 2020

For and on behalf of:

Shenzhen Precision Eyewear
Testing & Inspection Services Co., Ltd.

Manager: WenHua Li

Test Report

Report No.: PL2004308

Page 2 of 10

Tests Conducted

Note: The applicant's attention was drawn that the manufacturer should not use the frame materials which are known to cause irritation, allergic or toxic reaction during wear in a normal state of health against significant proportion of users.

Requirement			Testing		Results ¹	
Test Items		According to Clause		According to		
		EN	Clause	EN		Clause
General construction		166	6.1	--	--	P
Materials (Nickel release)		166	6.2	--	--	NA
Headbands		166	6.3	--	--	P
Field of vision		166	7.1.1	168	18	P
Refractive powers (Unmounted oculars covering one eye)	Spherical refractive powers	166	7.1.2.1.1	167	3.1	NA
	Astigmatic refractive powers					
	Prismatic refractive powers					
Refractive powers (Mounted oculars and covering both eyes)	Spherical refractive powers	166	7.1.2.1.2	167	3.2	Optical Class 1
	Astigmatic refractive powers					
	Prismatic refractive powers					
Transmittance	Oculars without filtering action	166	7.1.2.2.1	167	6	P
	Oculars with filtering action	166	7.1.2.2.2	167	6	NA
	Ultraviolet Filter	170	4	167	6	NA
	Sunglare Filter for Industrial Use	172	4.1	167	6	NA
Variations in transmittance (Exempt oculars without filtering action)	Oculars without corrective effect	166	7.1.2.2.3.1	167	7	NA
	Oculars with corrective effect	166	7.1.2.2.3.2	167	7	NA



Test Report

Report No.: PL2004308

Page 3 of 10

Tests Conducted

Requirement			Testing		Results		
Test Items		According to Clause		According to			
		EN	Clause	EN		Clause	
Diffusion of light			166	7.1.2.3	167	4	P
Quality of material and surface			166	7.1.3	167	5	P
Minimum robustness ²			166	7.1.4.1	168	4	NA
Increased robustness	Unmounted oculars		166	7.1.4.2.1	168	3.1	NA
	Complete eye-protectors and frame		166	7.1.4.2.2	168	3.2	P
Stability at an elevated temperature			166	7.1.5.1	168	5	P
Resistance to ultraviolet radiation (oculars only)			166	7.1.5.2	168	6	P
Resistance to corrosion (All metal parts only)			166	7.1.6	168	8	NA
Resistance to ignition			166	7.1.7	168	7	P
Protection against high-speed particles			166	7.2.2	168	9	NA
Protection against droplets and splashes			166	7.2.4	168	12	NA
Lateral protection			166	7.2.8	168	19	P
Resistance to fogging of oculars			166	7.3.2	168	16	P
Information supplied by the manufacturer			166	10	--	--	P

Remarks: 1. P = Pass; F = Fail; NA = Not Applicable; NR= Not require; X=Checked.

2. This requirement relates only to cover plates and oculars with filtering effect and not be assessed if these Items are intended to meet the requirements for increased robustness or resistance to high speed particles, in which case the requirements of 7.1.4.2 or 7.2.2 shall be met.

Test Report

Report No.: PL2004308

Page 4 of 10

Test Result

General construction — Clause 6.1/ Headbands – Clause 6.3

Sample No.	General construction		Headbands				Comment	Result(s)
	Defects		Width		Adjustable			
	Observed	Absent	Pass	Fail	Yes	No		
2004308-(01~03)		X	X		X		--	P
Requirements:								
1. General construction: Eye-Protectors shall be free of projections, sharp edges or other defects which are likely to cause discomfort or injury during use.								
2. Headbands: shall be at least 10mm wide over any portion which may come into contact with the wearer's head, and shall be adjustable or self-adjusting:								

Field of vision — Clause 7.1.1 / EN 168:2001 Clause 18

Sample No.	Head-form		Exhibit minimum field of vision defined in the standard		Comment	Result(s)
	Medium	Small	Yes	No		
2004308-(01~03)	X		X		--	P
Requirements:						
Eye-Protectors shall be exhibit field of vision an area of not less than 22 mm in the horizontal length and 20mm in the vertical width in front of each eye.						

Refractive powers— Clause 7.1.2.1 .2 / EN 167:2001 Clause 3.2

Sample No.	Refractive powers				Difference in prismatic refractive powers (cm/m)			Result(s)
	Spherical (m ⁻¹)		Astigmatic (m ⁻¹)		Horizontal		Vertical	
	Left	Right	Left	Right	Base Out	Base In		
2004308-01	+0.02	-0.02	0.00	0.04	0.00	--	0.01	Optical class 1
2004308-02	+0.01	-0.01	0.00	0.04	--	0.02	0.00	
2004308-03	+0.01	-0.01	0.00	0.00	0.00	--	0.00	
Requirement: Permissible tolerances for refractive powers :								
Optical class 1	±0.06		0.06		0.75	0.25	0.25	
Optical class 2	±0.12		0.12		1.00	0.25	0.25	
Optical class 3	+0.12~ -0.25		0.25		1.00	0.25	0.25	

Measurement Uncertainty (if necessary):

Test Report

Report No.: PL2004308

Page 5 of 10

Test Result

Transmittance (without filtering action) — Clause 7.1.2.2/ EN 167:2001 Clause 6

Sample No.	Requirements	Luminous Transmittance, T_v (%)		Result(s)
		Left	Right	
2004308-01	$T_v > 74.4\%$	91.1	92.2	P
2004308-02		90.3	91.6	P
2004308-03		90.2	91.5	P

Measurement Uncertainty (if necessary):

Diffusion of light — Clause 7.1.2.3 / EN 167:2001 Clause 4

Sample No.	Samples type	Diffusion of light (cd/m^2) / lx		Result(s)
		Left	Right	
2004308-01	III	0.30	0.20	P
2004308-02		0.16	0.07	P
2004308-03		0.12	0.08	P

Requirements:

The maximum value of the reduced luminance factor shall be :

- I 1.00(cd/m^2) / lx for welding filter;
- II 0.75(cd/m^2) / lx for oculars used in eye-protectors against high speed particles;
- III 0.50 (cd/m^2) / lx for all other oculars;

Measurement Uncertainty (if necessary):

Test Report

Report No.: PL2004308

Page 6 of 10

Test Result

Quality of material and surface — Clause 7.1.3 / EN 167:2001 Clause 5

Sample No.	Defects		Comment	Result(s)
	Observed	Absent		
2004308-(01~03)		X	--	P
Requirements: Except in a marginal area 5 mm wide, oculars shall be free from any significant defects likely to impair vision in use, such as bubbles, scratches, inclusions, dull spots, pitting, mould marks, scouring, grains, pocking, scaling and undulation.				

Increased robustness — Clause 7.1.4.2 / EN 168:2001 Clause 3.2

Sample No.	Test temperature (℃)	Test position	Defects		Comment	Result(s)
			Observed	Absent		
2004308-(04~05)	55	Left Frontal		X	--	P
2004308-06		Left Lateral		X	--	P
2004308-(07~08)	-5	Left Frontal		X	--	P
2004308-09		Left Lateral		X	--	P
2004308-(04~05)	55	Right Frontal		X	--	P
2004308-06		Right Lateral		X	--	P
2004308-(07~08)	-5	Right Frontal		X	--	P
2004308-09		Right Lateral		X	--	P
Requirements: The following defects shall not occur: 1. Ocular fracture: 2. Ocular deformation;						

Stability at an elevated temperature — Clause 7.1.5.1 / EN 168:2001 Clause 5

Sample No.	Apparent deformation		Comment	Result(s)
	Observed	Absent		
2004308-(04~06)		X	--	P
Requirements: Assembled eye-protectors shall show no apparent deformation;				

Test Report

Report No.: PL2004308

Page 7 of 10

Test Results

Resistance to ultraviolet radiation (oculars only) — Clause 7.1.5.2 / EN 168:2001 Clause 6

Samples type		Sample No.																				
Other oculars		2004308-01		2004308-02		2004308-03																
Test Items		Left	Right	Left	Right	Left	Right															
The relative change of luminous transmittance (%)	Before Expose	91.1	92.2	90.3	91.6	90.2	91.5															
	After Expose	90.8	92.0	89.0	91.6	90.7	91.4															
	Difference	-0.3	-0.2	-1.4	0.0	0.6	-0.1															
Reduced scattered light coefficient (cd/m²) / lx	Before Expose	0.30	0.20	0.16	0.07	0.12	0.08															
	After Expose	0.21	0.16	0.13	0.09	0.11	0.10															
Result(s)		P		P		P																
Requirements:																						
1. The relative change of luminous transmittance				2. Reduced scattered light coefficient																		
<table><tr><th colspan="2">Luminous transmittance</th><th rowspan="2">Permissible relative Change (%)</th></tr><tr><th>Less than (%)</th><th>Up to (%)</th></tr><tr><td>100</td><td>17.8</td><td>±5</td></tr><tr><td>17.8</td><td>0.44</td><td>±10</td></tr></table>			Luminous transmittance		Permissible relative Change (%)	Less than (%)	Up to (%)	100	17.8	±5	17.8	0.44	±10	<table><tr><td>The maximum value of the reduced luminance factor shall be:</td></tr><tr><td>- 1.00(cd/m²) / lx for welding filter;</td></tr><tr><td>- 0.75(cd/m²) / lx for oculars used in eye-protectors against high speed particles;</td></tr><tr><td>- 0.50 (cd/m²) / lx for all other oculars;</td></tr></table>					The maximum value of the reduced luminance factor shall be:	- 1.00(cd/m²) / lx for welding filter;	- 0.75(cd/m²) / lx for oculars used in eye-protectors against high speed particles;	- 0.50 (cd/m²) / lx for all other oculars;
Luminous transmittance		Permissible relative Change (%)																				
Less than (%)	Up to (%)																					
100	17.8	±5																				
17.8	0.44	±10																				
The maximum value of the reduced luminance factor shall be:																						
- 1.00(cd/m²) / lx for welding filter;																						
- 0.75(cd/m²) / lx for oculars used in eye-protectors against high speed particles;																						
- 0.50 (cd/m²) / lx for all other oculars;																						

Measurement Uncertainty (if necessary):

Resistance to ignition — Clause 7.1.7 / EN 168:2001 Clause 7

Sample No.	Continued combustion		Comment	Result(s)
	Yes	No		
2004308-(04~06)		X	--	P
Requirements:				
Eye-protectors shall be considered to be satisfactory if no parts ignites or continues to glow after removal of the steel rod.				

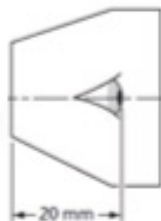
Test Report

Report No.: PL2004308

Page 8 of 10

Test Results

Lateral Protection — Clause 7.2.8 / EN 168:2001 Clause 19

Sample No.	Lateral region coverage		Comment	Result(s)
	Pass	Fail		
2004308-(01~03)	X		---	P
Requirements: Eye-protectors claimed to lateral protection shall pass the lateral region coverage. These region being a 20mm wide with 10 mm radial ends struck from the front and side impact points (See figure 1)			Figure 1: 	

Resistant to Fogging of Oculars — Clause 7.3.2/ EN 168:2001 Clause 16

Sample No.	Position	After Fogging Test Evaluation		Comment	Result
		Fogging	Free Fogging		
2004308-10	Left		X	--	P
	Right		X	--	P
Requirements: If oculars are described as resistant to fogging they shall remain free from fogging for a minimum of 8s.					



Test Report

Report No.: PL2004308

Page 9 of 10

Test Results

Information supplied by the manufacturer — Clause 10

The manufacturer shall be provided with each eye-protector at least the following information:

Information	Results
Name and address of the manufacturer	P
The number of this standards	P
The eye-protector model identification	P
Instructions for storage, use and maintenance	P
Specific instructions for cleaning and disinfection	P
Details of the filed of use, protection capabilities and performance characteristics	P
The significance of making on the fame and the ocular	P
A warning concerning the compatibility of marking	P
A warning that materials which may come into contact with the wearer's skin could cause allergic reactions to susceptible individuals	P
A warning that scratched or damaged oculars should be replaced	P
A warning that eye-protection against high speed worn over standard ophthalmic spectacles may transmit impacts, thus creating a hazard to the wearer.	P

Test Report

Report No.: PL2004308

Page 10 of 10

Sample Photo:



-----Report End-----

The test report is effective only with both signature and specialized stamp. The result(s) shown in this report refer only to the sample(s) tested. Without written approval of PEL (Precision Eyewear Testing & Inspection Co., Ltd.), this report can't be reproduced except in full.

