

**HEAD OFFICE**

Owen Mumford Ltd  
Brook Hill, Woodstock  
Oxfordshire OX20 1TU, United Kingdom  
**T** : +44 (0)1993 812021  
**E** : info@owenmumford.co.uk

**UK / IRELAND SALES & MARKETING**

Owen Mumford Ltd  
Primsdown Industrial Estate  
Worcester Road, Chipping Norton  
Oxfordshire OX7 5XP, United Kingdom  
**T** : +44 (0)1608 645555  
**E** : info@owenmumford.co.uk

**NORTH AMERICA**

**T** : +1 770-977-2226  
**E** : info@owenmumfordinc.com

**FRANCE**

**T** : +33 (0)2 32 51 88 70  
**E** : fr.info@owenmumford.com

**GERMANY**

**T** : +49 (0)6026/97750  
**E** : de.info@owenmumford.com

**CHINA**

**T** : +86 (0)21 6032 3528  
**E** : cn.info@owenmumford.com

**MALAYSIA**

**T** : +60 (0)7 553 9617  
**E** : my.info@owenmumford.com

**MEXICO**

**T** : +52 (55) 3733 8759  
**E** : mx.info@owenmumford.com

[owenmumford.com](http://owenmumford.com)

**OUR ACCREDITATIONS INCLUDE:**

- ISO13485 Certification
- Medical Device Directive 93/42/EEC
- CMDCAS 13485 (Canadian Health Authority)
- Environmental Management Systems ISO14001



© Copyright Owen Mumford 2016  
Ufouch-INTL/IMCC/0116/02

 **ntuitive**

The simple, easy safety lancet



**Unistik® Touch**

# Comfort

## Comfort Zone Technology®

Unistik® Touch features our patented Comfort Zone Technology®. This simply means that when you press the lancet against the skin, the 8 raised dots masks the pain sensation created by the needle.<sup>1</sup>

Comfort Zone Technology® is based on the 'Gate Theory'<sup>1</sup> and works by sending a 'comfort' message to the brain which interferes with the pain message produced during the lancing.

Simply put: Unistik® means comfort first.



## Clinical Evidence<sup>2</sup>

Lancet trials showed that Comfort Zone Technology® (CZT®) reduced the overall pain experienced during the blood sampling procedure when compared to devices without Comfort Zone Technology®.


1. Melzack, R. and Wall, P.D. (1982) *The challenge of pain*. Harmondsworth, UK: Penguin, p233. 2. Dewland, P. and Edwards, C. (2007) *A single-blind, randomised, 8-way crossover study to compare the blood volume and pain perception of capillary blood sampling*, UK.

# Simple

## 3 easy steps

Unistik® Touch is designed to make blood sampling simple and easy. Touch activation means the lancet can be used in 3 easy steps.


**1.**

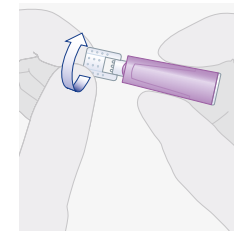

### TWIST

Twist off and discard lancet cap. **DO NOT PULL.**


**2.**

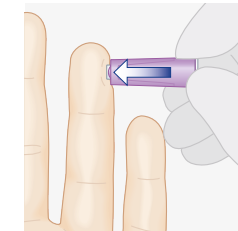

### CLICK

Press platform firmly against sample site. Lancet will automatically retract into body of device.

**3.**

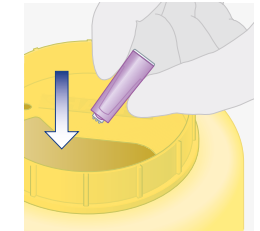

### DISPOSE

Dispose of in a suitable sharps container.

Unistik® Touch is available in four different needle sizes to allow you to choose a suitable lancet for every test you perform.



For frequent use and paediatrics



**30G**  
1.5mm

100: 000361  
200: 000362



Frequent blood glucose testing



**28G**  
1.8mm

100: 000363  
200: 000364



Blood glucose and standard testing



**23G**  
2.0mm

100: 000365  
200: 000366



Larger samples



**21G**  
2.0mm

100: 000367  
200: 000368